
Scheme 1

Solid Phase Protein Synthesis

Native Chemical Ligations in an N- to C-Terminal Direction

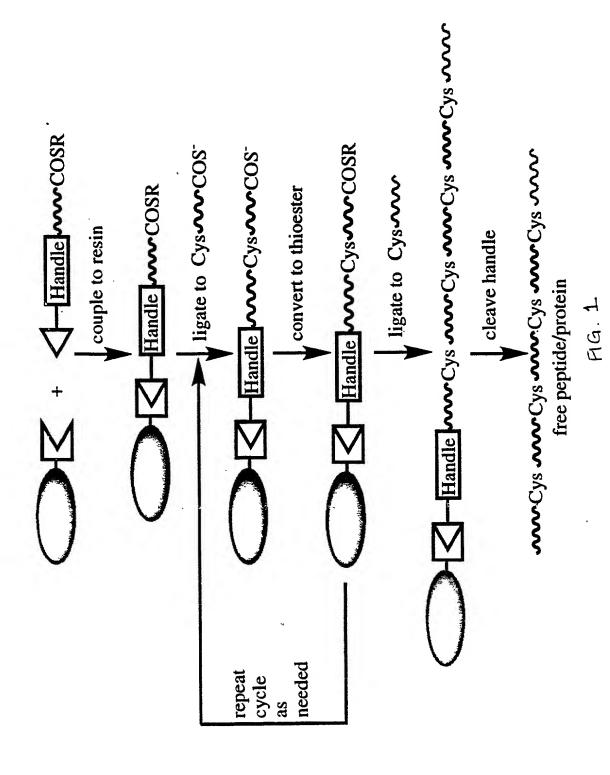

Cys + COSR Stability Under Ligation Conditions

In the absence of a thioester peptide
 H-CGFRVREFGDNTA-COSH MW=1487.6
 6M GU-HCL, 0.1M NaPi, 0.5% thiophenol, room temperature, overnight

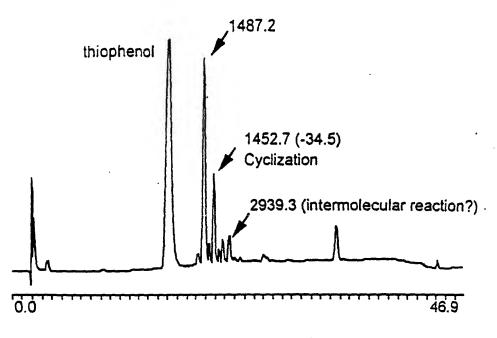


FIG. 2A

2. In the presence of a thioester peptide

H-CGFRVREFGDNTA-COSH MW=1487.6 + H-DSVISLSGDH-SPAL MW= 1230.2

MW of Ligation product = 2498.7

6M GU-HCL, 0.1M NaPi, 0.5% thiophenol, room temperature, overnight

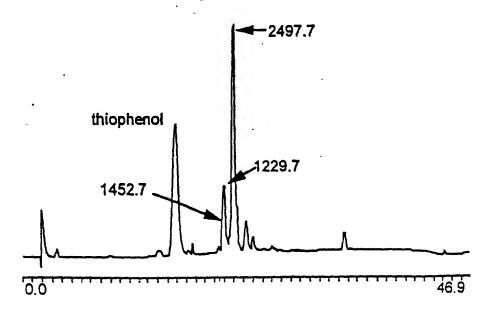



FIG. 2B

MSC Removal Experiments

MSC-CTSAGPHFNPLSRKHG-OH MW=1859.1

H-CTSAGPHFNPLSRKHG-OH MW=1708.9


Aliquot of peptide in 6M Gu•HCl, 0.1M NaPi, pH 7.5 was diluted into 1N NaOH for two minutes, quenched with 1N HCl

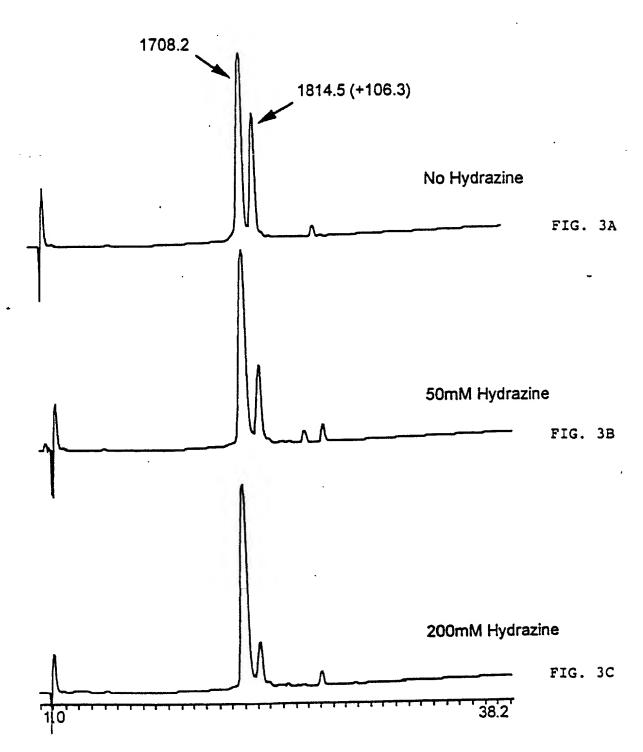

MSC Removal Experiments (Cont'd)


Lev-MSC-LTEGLHGFHVHEFGDNTAGCTSAGPHFNPLSRKHG-COSH MW=4022.4 H-LTEGLHGFHVHEFGDNTAGCTSAGPHFNPLSRKHG-COSH MW=3745.1

Aliquot of peptide in 6M Gu•HCl, 0.1M NaAc, pH 4.6 was diluted into 6M Gu•HCl, 0.1M NaAc, pH 14 for two minutes, quenched with 6M Gu•HCl, 0.1M NaAc, pH 2.0

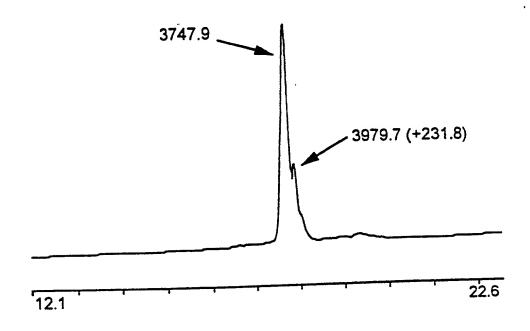

Resin Preparation

Resin Preparation

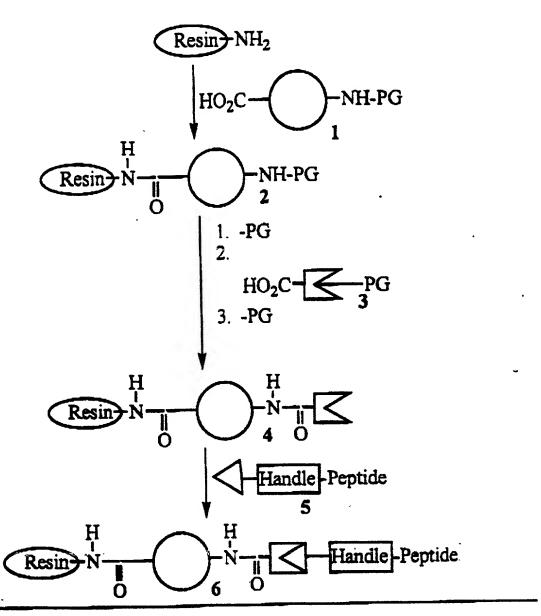

= cleavable linker used for monitoring with Maldi, Electrospray Mass Spec, etc...

PG = protecting group

HO₂C-

= functional group added to resin to couple with peptide

Handle-Peptide

= peptide functionalized with 1) a cleavable handle for release of peptide/protein from the resin at completion of synthesis and 2) functional group to couple to resin

Derivatization of Segment 1 (N-terminal)


Boc-LTEGLHGFHVHEFGDNTAGCTSAGPHFNPLSRKHG-COS(AM)

Poly mer-Supported Ligation on PEGA

Lev-MSC-LTEGLHGFHVHEFGDNTAGCTSAGPHFNPLSRKHG-COSH (1) + Resin-PCL-ONH2

 \downarrow 1. pH 4.6, 6M Gu-HCl, 0.1 acetate

Resin-PCL-oxime-MSC-LTEGLHGFHVHEFGDNTAGCTSAGPHFNPLSRKHG-COSH (1)

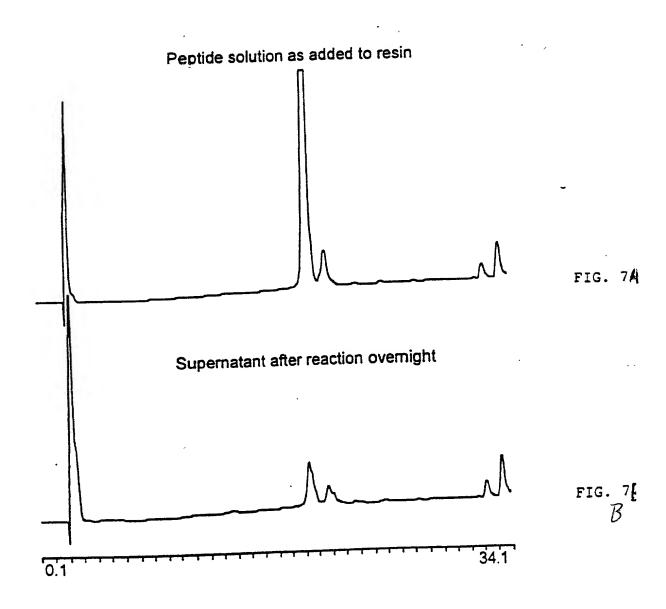


Polymer-Supported Ligation on Isto

Lev-MSC-LTEGLHGFHVHEFGDNTAGCTSAGPHFNPLSRKHG-COSH (1) + Resin-PCL-ONH2

↓ 1. pH 4.6, 6M Gu•HCl, 0.1 acetate

Resin-PCL-oxime-MSC-LTEGLHGFHVHEFGDNTAGCTSAGPHFNPLSRKHG-COSH (1)

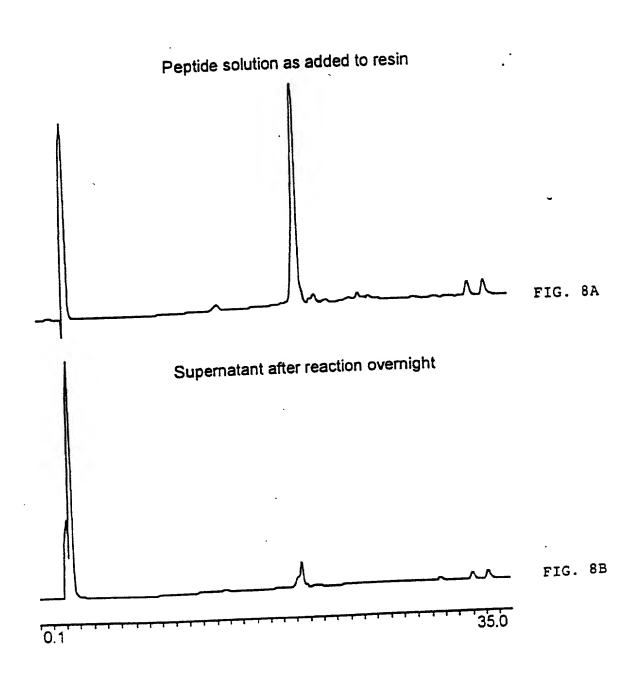


Polymer-Supported Ligation on Isco

Lev-MSC-LTEGL: GFHVHEFGDNTAGCTSAGPHFNPLSRKHG-COSH (1) + Resin-PCL-ONH2

↓ 1. pH 4.6, 6M Gu•HCl, 0.1 acetate

Resin-PCL-oxime-MSC-LTEGLHGFHVHEFGDNTAGCTSAGPHFNPLSRKHG-COSH (1)
Maldi Mass = 4022, Base Cleavage Mass = 3745

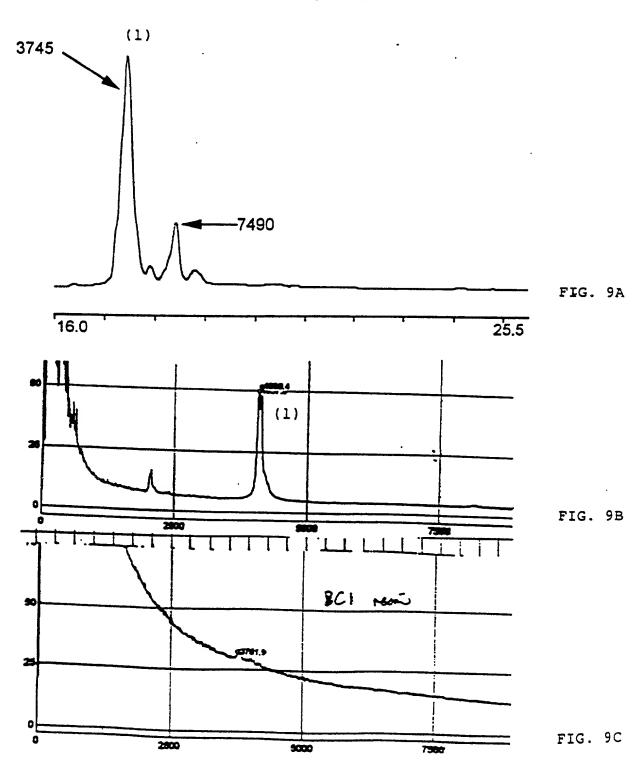

Polymer-Supported Ligation 67 Isco

Resin-PCL-oxime-MSC-Li EGLHGFHVHEFGDNTAGCTSAGPHFNPLSRKHG-COSAc (1)


Maldi Mass = 4080, Base Cleavage Mass = 3729


+ H-CGFRVREFGDNTA-COSH (2)

↓ 3. pH 7.5, 6M Gu•HCl, 0.1M phosphate, 0.5% thiophenol

Resin-PCL-oxime-MSC-LTEGLHGFHVHEFGDNTAGCTSAGPHFNPLSRKHGCGFRVREF-GDNTA-COSH (1+2)

Maldi Mass = 5476, Base Cleavage Mass = 5199

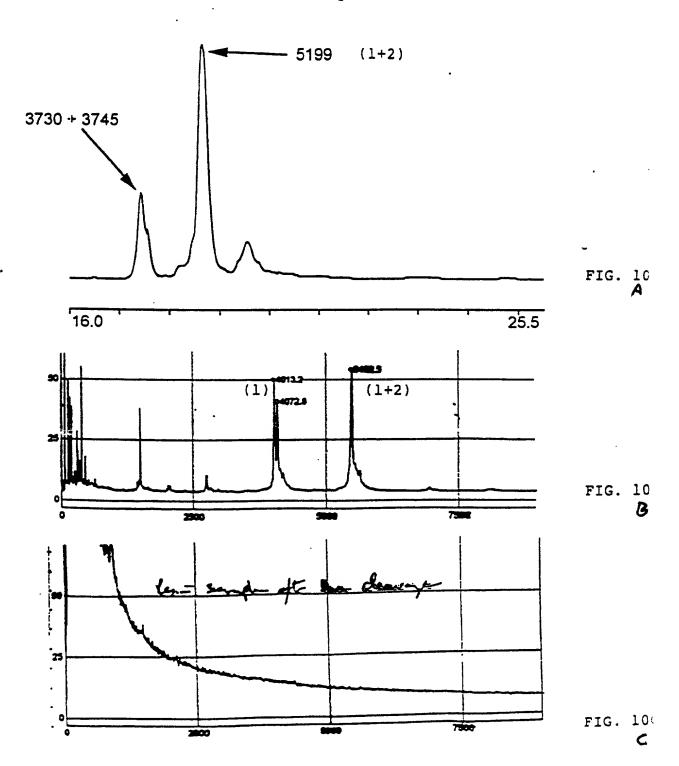

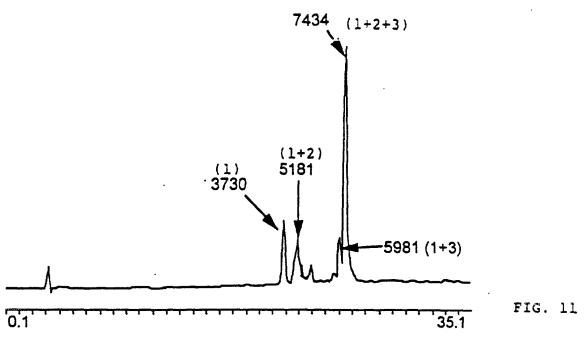

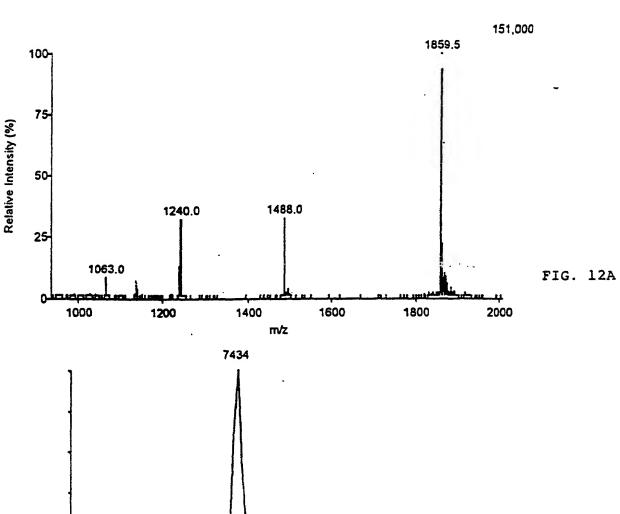
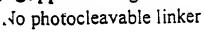


FIG. 12B

Though the time of the thing then the term of the term

i=1

4 W T W



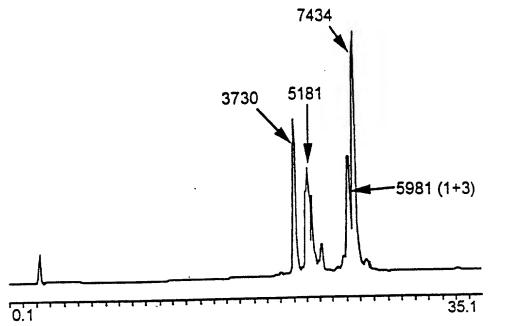


FIG. 13

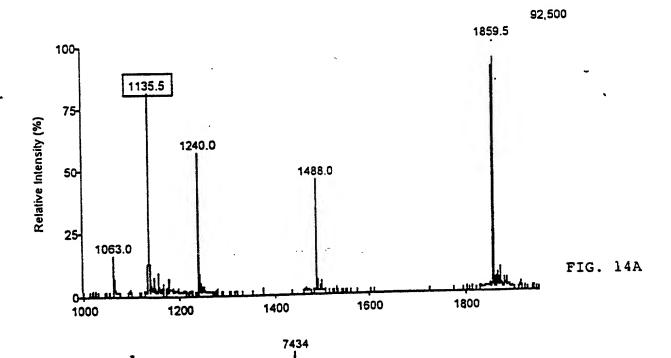
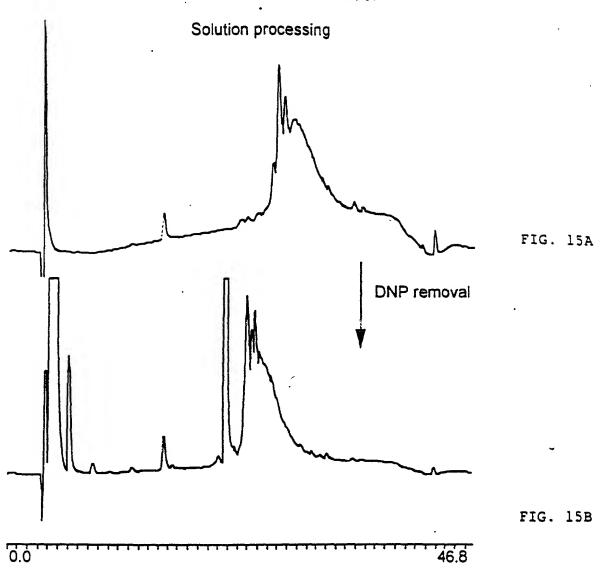


FIG. 14B

On Resin Purification



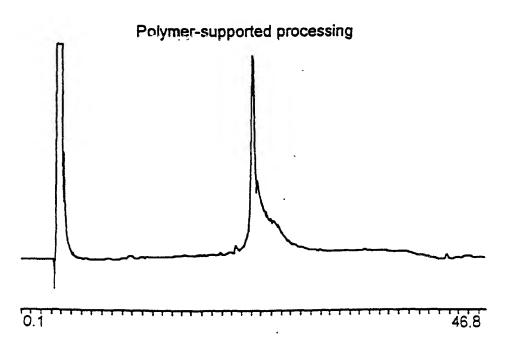


FIG. 15C

Synthesis of MIF by Solid Phase Native Ligations

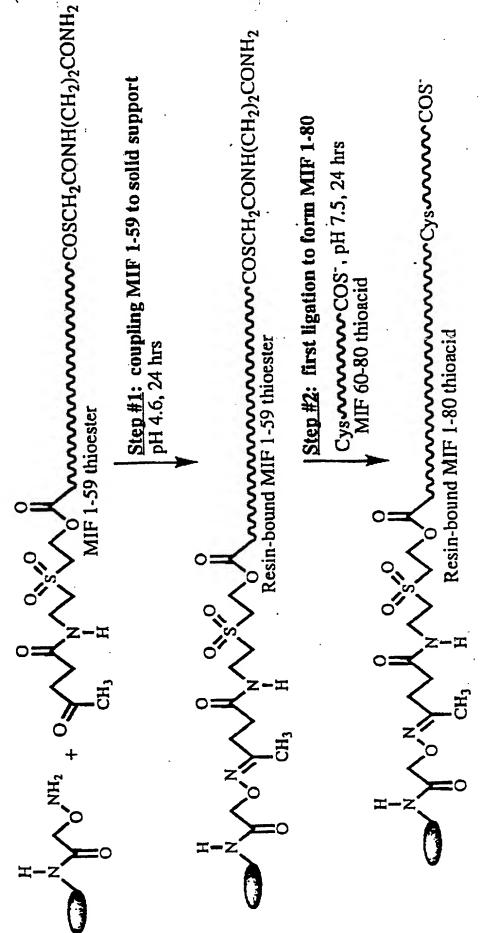


FIG. 164

Step #3: conversion to thioester BrCH₂CO₂H, pH 4.6, 15 min

www.cosch.co.h Resin-bound MIF 1-80 thioester

Step #4: second ligation to form MIF 1-115

Cys. MIF 81-115

Resin-bound MIF 1-115


Step #5: base cleavage from solid support pH ~14, 2 min

www.cyswww.cyswww.

Free MIF 1-115

FIG. 168

Modification of N-terminal Peptide Segment and Solid Support

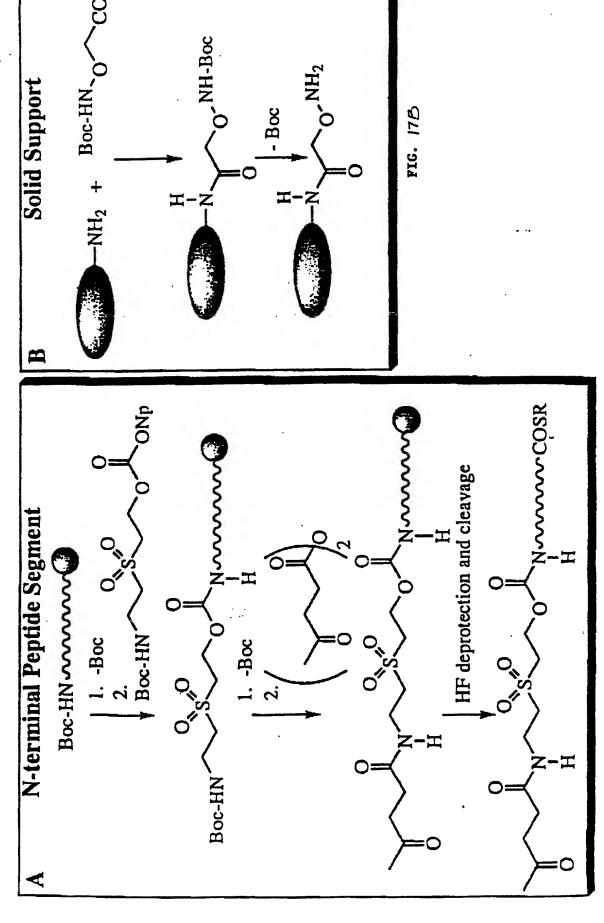


Coupling of MIF 1-59 to Solid Support

Ketone-MSC handle-Met ¹-MIF 2-58-Leu⁵⁹-SAc-βAla-CO₂H

#1

Isco Oxime-MSC handle-Met¹-MIF 2-58-Leu⁵⁹-SAc-βAla-CO₂H Expected base cleavage mass = 6271


FIG. (BA

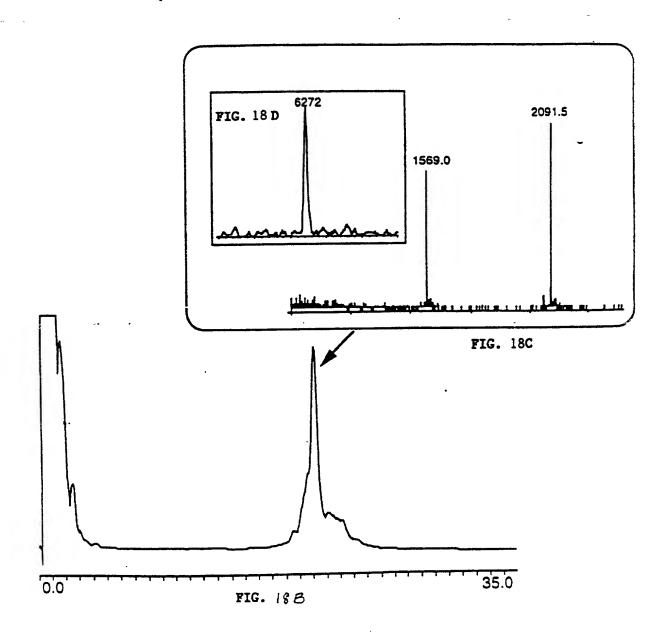

#2 Cys⁶⁰-MIF 61-79-Leu³⁰-COSH

Isco Oxime-MSC handle-Met¹-MIF 2-79-Leu⁸⁰-COSH Expected base cleavage mass = 8502

FIG. 19A

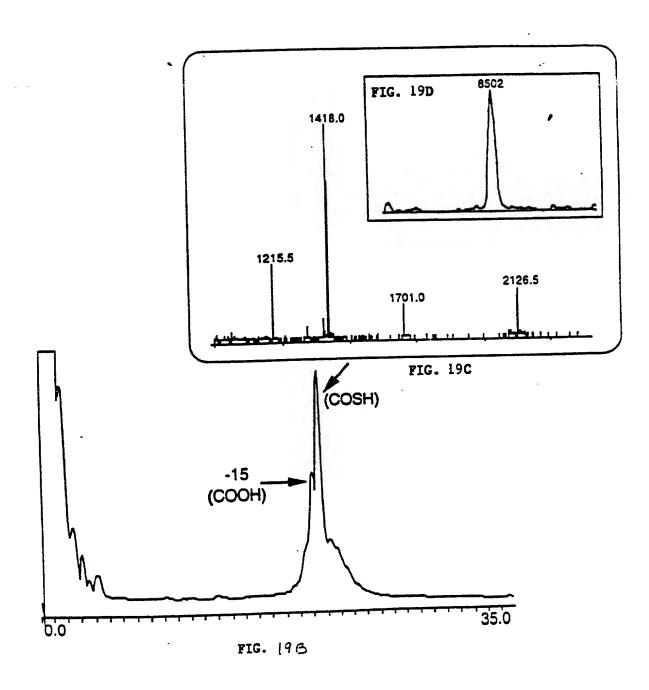


Isco Oxime-MSC handle-Met¹-MIF 2-79-Leu⁸⁰-COSAc


Cys⁸¹-MIF 82-114-Ala¹¹⁵-CO₂H
6M Gu•HCl, 0.1 M Na Pi, 0.5% thiophenol
0.15 M Methionine, pH 7.5

Isco — Oxime-MSC handle-Met¹-MIF 2-114-Ala¹¹⁵-CO₂H Expected base cleavage mass = 12450

FIG. ZOA

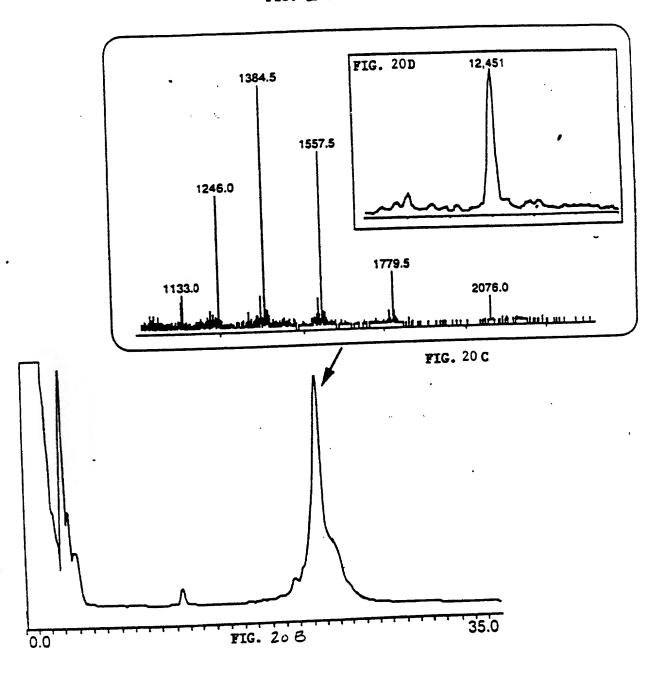

Solid Phase Chemical Ligations in the C- to N-terminal Direction

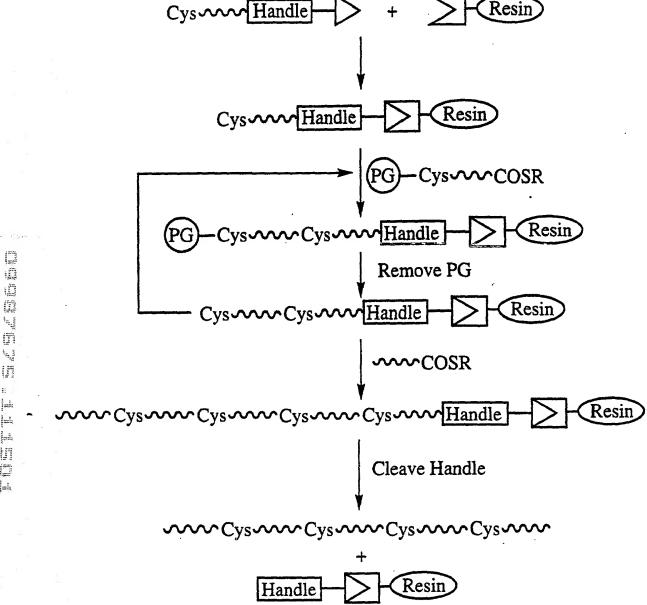

Solid Phase Chemical Ligations in the C- to N-Terminal Direction Synthesis of Phospholipase A2, Group 5 (PLA2G5)


Figure 22

Synthesis of Cam ester derivative


Synthesis of C-terminal peptide segment

1 1877

Universal Solid Phase Chemical Ligation (Bidirectional Ligations: C- to N-Terminal Ligations First)

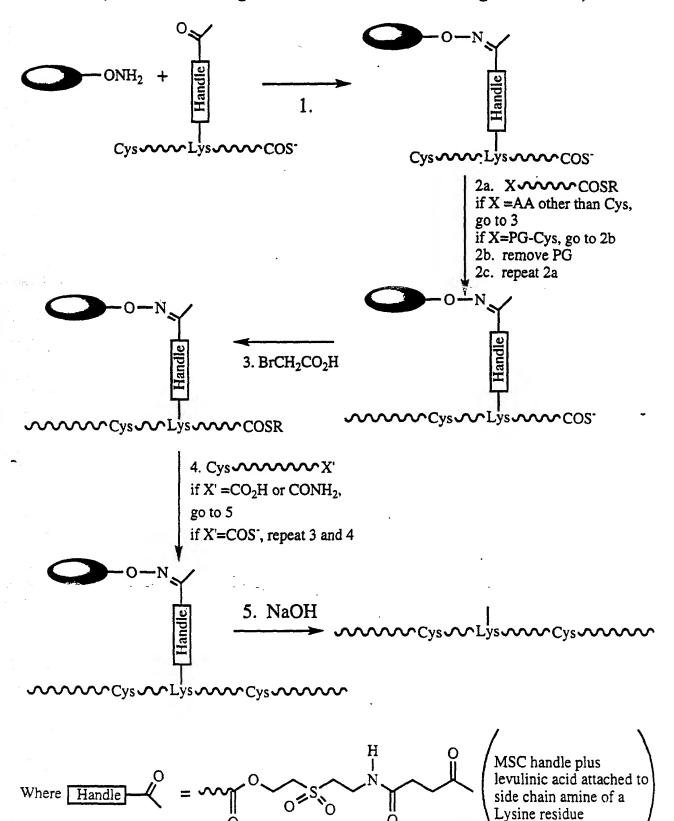


ĺ.

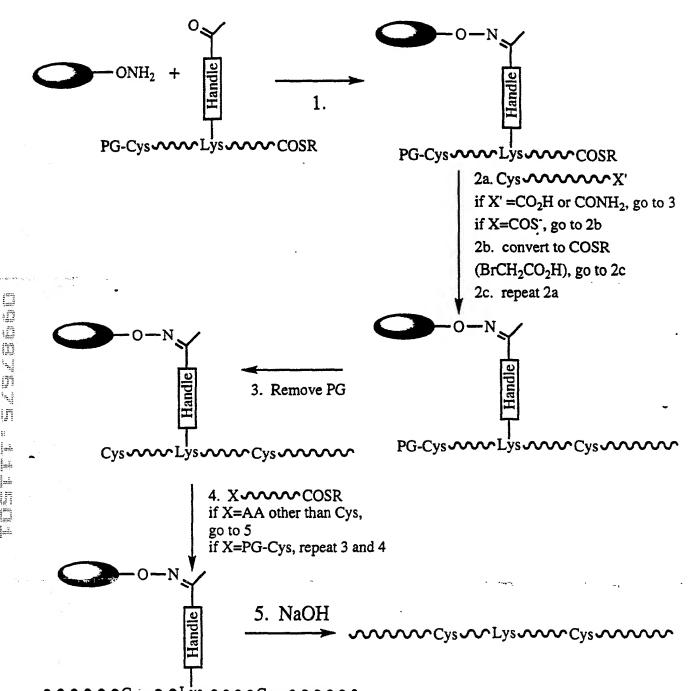


Synthesis of Modified Peptide Segment for Universal Solid Phase Chemical Ligation

Starting with an appropriate resin (thioester or thioacid generating), synthesize the peptide using standard Boc protocols until the Lys residue of choice is reached. Couple a Boc-Lys(Fmoc)-OH, then continue the rest of the synthesis.

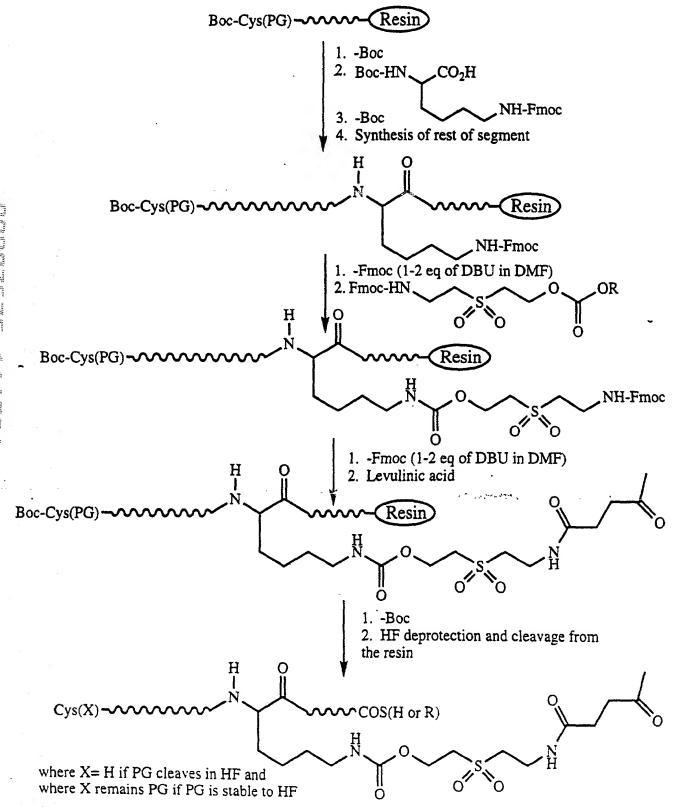

Synthesis of C5a by Solid Phase Chemical Ligations in the N- to C-Terminal Direction

1 21 47
TLQKKIEEIAAKYKHSVVKK**CCYDGACVNNDETCEQRAARISLGPK**CIKAFTECC
VVASQLRANISHKDMQLGR
74

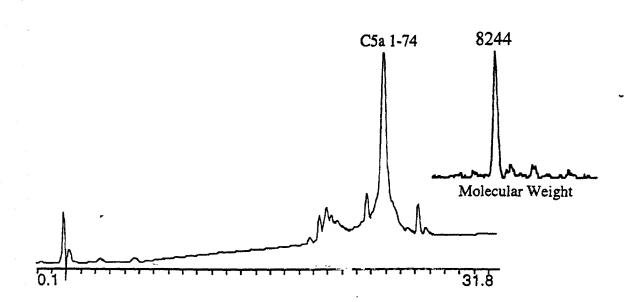


Synthesis of C-terminal peptide segment

F1G1. 27

Synthesis of a Random Sequence by Solid Phase Chemical Ligations in the C- to N-terminal Direction


ALTKYGFYGCYGRLEEKG**CADRKNILA** 19 27 10

Using Fmoc Protection

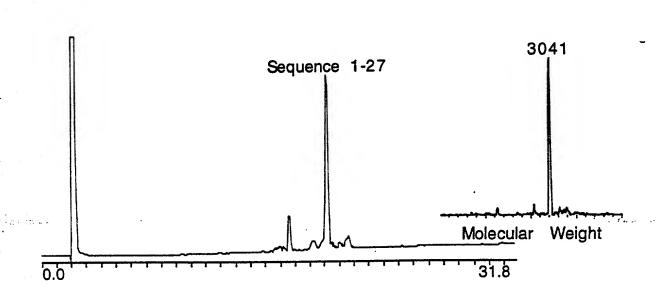


Synthesis of a Random Sequence by Solid Phase Chemical Ligations in the C- to N-terminal Direction Using ACM Protection

ALTKYGFYGCYGRLEEKGCADRKNILA
1 10 19 27

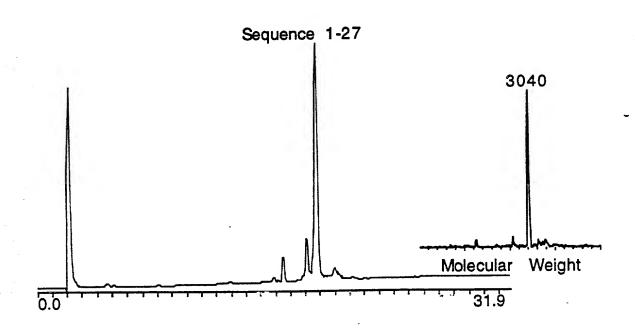

Synthesis of Phospholipase A2 Group 5 by Solid Phase Chemical Ligations in the C- to N-Terminal Direction

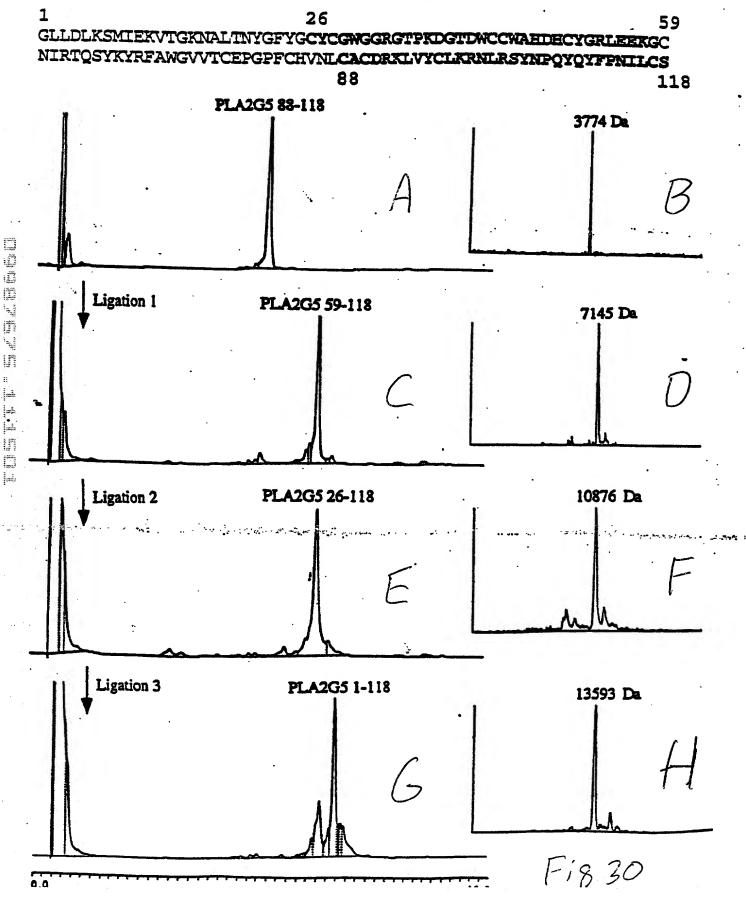

Synthesis of Phospholipase A2 Group 5 by Solid Phase Chemical Ligations in the C- to N-Terminal Direction

